

# AM4515ZT - EDGE

**Brand:** Dino-Lite  
**Product Code:** AM4515ZT - EDGE



## Short Description

1.3 Megapixel Edge sensor  
 Versatility due to exchangeable caps  
 20-220x magnification  
 Polarizer (anti-reflection)  
 Automatic Magnification Reading

## Description

The Dino-Lite Edge models provide superior image quality and a high level of flexibility due to the adaptable cap design. The built-in adjustable polarizer reduces glare and reflection on shiny objects.



USB 2.0



Measurement  
functionality



Adjustable  
~20x-220x



Automatic  
Magnification  
Reading (AMR)



Magnification  
lock




1.3 megapixel  
1280 x 1024



8 White LEDs  
switchable



Polarizer  
anti-reflection

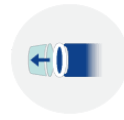

Exchangeable  
cap

With the built-in automatic magnification reading, measurements can be performed easily and quickly. Different caps are supplied with the microscope making the Dino-Lite Edge

the perfect choice for any professional application.

*Also available: AM4515T - Edge without Polarizer*

**Working distance/field of view/depth of field**

MAGNIFICATION RATE	WORKING DISTANCE 1	WORKING DISTANCE 2*	FIELD OF VIEW(X)	FIELD OF VIEW(Y)	DEPTH OF FIELD
<b>Listed values may differ slightly</b>		<b>* Without front cap</b>			<b>Unit = mm</b>
20	52.7	60.2	19.5	15.6	2.5
30	26.0	33.5	13.0	10.4	1.8
40	13.4	20.9	9.8	7.8	1.5
50	6.4	13.9	7.8	6.3	-
60	2.2	9.7	6.5	5.2	-
70	-	7.1	5.6	4.5	1.0
80	-	5.5	4.9	3.9	-
90	-	4.5	4.3	3.5	-
100	-	4.1	3.9	3.1	-
110	-	4.0	3.6	2.8	-
120	-	4.1	3.3	2.6	-
130	-	4.5	3.0	2.4	-
140	-	5.0	2.8	2.2	-
150	-	5.6	2.6	2.1	-
160	-	6.3	2.4	2.0	-
170	-	7.1	2.3	1.8	-
180	0.5	8.0	2.2	1.7	-
190	1.4	8.9	2.1	1.6	-
200	2.4	9.9	2.0	1.6	-
210	3.4	10.9	1.9	1.5	-
220	4.4	11.9	1.8	1.4	0.1

**Specification**

Lighting	
Light/ LED type	White

Number of LEDs	8
LED on/off switchable:	Yes
Infrared filter	IR cut-filter >650 nm
Diffuser available	Yes (N3C-D included)
Emission filter	No
Polarizer	Yes, linear
<b>Optics</b>	
Magnification	20-220x
Macro zoom	No
Working distance	Standard
Lens type	Glass with anti-reflection coating
<b>Sensor</b>	
Sensor type	CMOS
Resolution	1.3 Megapixel (1280x1024)
Maximum frame rate	30 fps
<b>Compatibility</b>	
Interface	USB 2.0
Operating system	Windows 7,8,10 & 11, MacOS 10.9 and up
Included software	DinoCapture 2.0 (Windows), DinoXscope (Mac OS)
Supported image formats (Windows)	BMP, GIF, PNG, JPG, TIF, RAS, PNM, TGA, PCX, MNG, WBMP, JP2, JPC, PGX
Supported video formats (Windows)	WMV, FLV, SWF
Imaging standards	DirectShow, UVC
Wifi	Wireless-ready, requires the WF-10 WiFi streamer (optional)
<b>Housing</b>	
Housing material	Composite/ plastic housing
Magnification lock	Yes
Dimensions	10.5cm (L) x 3.2cm (D)
Weight	105g
Cable length	1.8m
<b>Features</b>	
Special feature	Automatic Magnification Reading (AMR)
Measurement	Yes
Calibration	Yes
Microtouch sensor	Yes
ESD safe	No
<b>Information</b>	
Package contents	Microscope, carry pouch, software CD, calibration target, user manual, N3C-O- Open cap, N3C-C-

	Closed cap, N3C-D- diffuser cap, N3C-E- extension cap, N3C-L- Long cap, N3C-S- Side light cap
Warranty information	2 years European warranty
Regulatory approval	CE, FCC, ROHS
Price range	€600,00 - €800,00

## Product Gallery

