

AM7115MZT - EDGE

Brand: Dino-Lite
Product Code: AM7115MZT - EDGE

Short Description

5 Megapixel resolution
 Robust metal housing
 20-220x magnification
 Polarizer
 Flexible LED Control (FLC)

Description

With the use of the latest, cutting-edge optics, a brand new 5 megapixel sensor and several special features such as Flexible LED Control (FLC) and full aluminum body, the Dino-Lite AM7115MZT is a marvel of technology and the best choice for the high-demanding professional. The Dino-Lite AM7115MZT offers superb image quality and color reproduction in a robust, compact and appealing housing.



USB 2.0



Measurement
functionality



Adjustable
~20x-220x



Magnification
lock



Metal
housing



5 megapixel
2592 x 1944



8 White LEDs
switchable




Standard working
distance



Polarizer
anti-reflection



Flexible LED
Control FLC

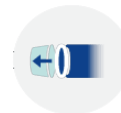

Exchangeable
cap

The AM7115MZT is part of the universal range and with its standard working distance and wide magnification range from 20 – 220x, it is suitable for a wide range of applications. Because of the built-in polarization filter this model is ideal when working with shiny or reflective objects such as metal, plastic, glass, jewelry, electronics, etc.

The main features of the AM7115MZT are:

- 5 Megapixel Edge sensor
- 20-220x magnification
- Integrated polarizer
- Aluminum housing
- Flexible LED Control (FLC)
- Extensive measurement functions
- Calibration
- Exchangeable front caps
- And more...

Working distance/field of view/depth of field

MAGNIFICATION RATE	WORKING DISTANCE 1	WORKING DISTANCE 2*	FIELD OF VIEW(X)	FIELD OF VIEW(Y)	DEPTH OF FIELD
Listed values may differ slightly		* Without front cap			Unit = mm
20	52.7	60.2	19.5	15.6	2.5
30	26.0	33.5	13.0	10.4	1.8
40	13.4	20.9	9.8	7.8	1.5
50	6.4	13.9	7.8	6.3	-
60	2.2	9.7	6.5	5.2	-
70	-	7.1	5.6	4.5	1.0
80	-	5.5	4.9	3.9	-
90	-	4.5	4.3	3.5	-
100	-	4.1	3.9	3.1	-
110	-	4.0	3.6	2.8	-
120	-	4.1	3.3	2.6	-
130	-	4.5	3.0	2.4	-
140	-	5.0	2.8	2.2	-
150	-	5.6	2.6	2.1	-
160	-	6.3	2.4	2.0	-
170	-	7.1	2.3	1.8	-
180	-	8.0	2.2	1.7	-
190	-	8.9	2.1	1.6	-
200	-	9.9	2.0	1.6	-
210	-	10.9	1.9	1.5	-

Specification

Lighting	
Light/ LED type	White
Number of LEDs	8
LED on/off switchable:	Yes
Infrared filter	IR cut-filter >650 nm
Diffuser available	Yes (N3C-D included)
Emission filter	No
Polarizer	Yes, linear
Optics	
Magnification	20-220x
Macro zoom	No
Working distance	Standard
Lens type	Glass with anti-reflection coating
Sensor	
Sensor type	CMOS
Resolution	5 Megapixel (2592x1944)
Maximum frame rate	30 fps
Compatibility	
Interface	USB 2.0
Operating system	Windows 7, 8, 10 & 11, MacOS 10.9 and up
Included software	DinoCapture 2.0 (Windows), DinoXcope (Mac OS)
Supported image formats (Windows)	BMP, GIF, PNG, JPG, TIF, RAS, PNM, TGA, PCX, MNG, WBMP, JP2, JPC, PGX
Supported video formats (Windows)	WMV, FLV, SWF
Supported image formats (MacOS)	JPEG, PNG
Supported video formats (MacOS)	MOV (max 1.3MP)
Imaging standards	DirectShow, UVC
Wifi	Wireless-ready, requires the WF-10 WiFi streamer (optional)
Housing	
Housing material	Metal housing
Dimensions	10.5cm (L) x 3.2cm (D)
Weight	135gr
Cable length	1.8m
Features	

Special feature	Flexible LED Control (FLC)
Measurement	Yes
Calibration	Yes
Microtouch sensor	Yes
ESD safe	Yes
Information	
Package contents	Microscope, carry pouch, software CD, calibration target, user manual, N3C-O- Open cap, N3C-C- Closed cap, N3C-D- diffuser cap, N3c-extension cap, N3C-L- Long cap, N3C-S- Side light cap
Warranty information	2 years European warranty
Regulatory approval	CE, FCC, ROHS
Price range	€750,00 - €1000,00

Product Gallery



