

AM4113ZT4

Brand: Dino-Lite
Product Code: AM4113ZT4



Short Description

1.3 Megapixel resolution
With magnification lock
400-470x magnification
Polarizer (anti-reflection)

Description

The AM4113ZT4 with a high magnification of 400~470x reveals even the smallest details on all kinds of objects. The built-in, switchable polarizer makes the AM4113ZT4 the right choice for working with high glare materials such as jewelry, plastics and metals.



USB 2.0



Measurement
functionality



Adjustable
~400x -470x



Magnification
lock




1,3 Megapixels
1280 x 1024



8 White LEDs
switchable



Standard working
distance

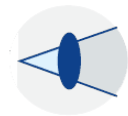

Polarizer
anti-reflection

Working distance/field of view/depth of field

MAGNIFICATION	WORKING DISTANCE	FIELD OF VIEW(X)	FIELD OF VIEW(Y)	DEPTH OF FIELD
430	-0.45	0.91	0.73	0.1
440	0.4	0.89	0.71	0.1
450	1.35	0.87	0.7	0.1
460	2.3	0.85	0.68	0.1
470	3.15	0.83	0.67	0.1

Listed values may differ slightly

Unit = mm

Specification

Lighting	
Light/ LED type	White
Number of LEDs	8
LED on/off switchable:	Yes
Infrared filter	IR cut-filter >650 nm
Diffuser available	No
Emission filter	No
Polarizer	Yes, linear
Optics	
Magnification	400-470x
Macro zoom	No
Working distance	Standard
Lens type	Glass with anti-reflection coating
Sensor	
Sensor type	CMOS, UKCA – available upon request – marketing@dino-lite.eu
Resolution	1.3 Megapixel (1280x1024)
Maximum frame rate	30 fps
Compatibility	
Interface	USB 2.0
Operating system	Windows 7,8,10 & 11, MacOS 10.9 and up
Included software	DinoCapture 2.0 (Windows), DinoXcope (Mac OS)
Supported image formats (Windows)	BMP, GIF, PNG, JPG, TIF, RAS, PNM, TGA, PCX, MNG, WBMP, JP2, JPC, PGX
Supported video formats (Windows)	WMV, FLV, SWF
Supported image formats (MacOS)	JPEG, PNG
Supported video formats	MOV

(MacOS)	
Imaging standards	DirectShow, UVC
Wifi	Wireless-ready, requires the WF-10 WiFi streamer (optional)
Housing	
Housing material	Composite/ plastic housing
Magnification lock	Yes
Dimensions	10.5cm (L) x 3.2cm (D)
Weight	105g
Cable length	1.8m
Features	
Special feature	No
Measurement	Yes
Calibration	Yes
Microtouch sensor	Yes
ESD safe	No
Information	
Package contents	Microscope, carry pouch, software CD, calibration target, ZT-Z-CC1 cap closure, user manual
Warranty information	2 years European warranty
Regulatory approval	CE, FCC, ROHS, UKCA – available upon request – marketing@dino-lite.eu
Price range	€400,00 - €600,00