

AM73515MZT

Brand: Dino-Lite
Product Code: AM73515MZT

Short Description

5 Megapixel Edge sensor
USB 3.0, maximum 45fps
10-220x magnification
Automatic Magnification Reading (AMR)
Flexible LED Control (FLC)
Polarizer

Description

With the use of the USB 3.0 connection the AM73515MZT is able to capture smooth and high-quality images. The Edge series optics and the special feature Flexible LED Control (FLC) and Automatic Magnification Reading (AMR) makes the AM73515MZT the best choice for the more demanding applications.

USB 3.0 adds a new transfer mode called "SuperSpeed" (SS) capable of transferring data up to 5Gbits/s (625MB/s). Dino-Lite models with USB 3.0 offer an image transfer speed of up to 45FPS at a 1280 x 960 resolution. SuperSpeed allows you to take advantage of the full power of the Dino-Lite Edge Series optics including improved colour accuracy and image quality. The AMR function automatically detects and displays the magnification of the Dino-Lite. The magnification is displayed within the Dino-Lite software and stored with the captured picture and simplifies the measuring process. The AM73515MZT is one of the most high-end Dino-Lite products available and the top-level product within the Dino-Lite USB 3.0 range. For professionals who are looking for high speed data transfer, the best available image quality and are looking for advanced features this is the model of choice.




USB 3.0



High speed



Adjustable
~10-220x

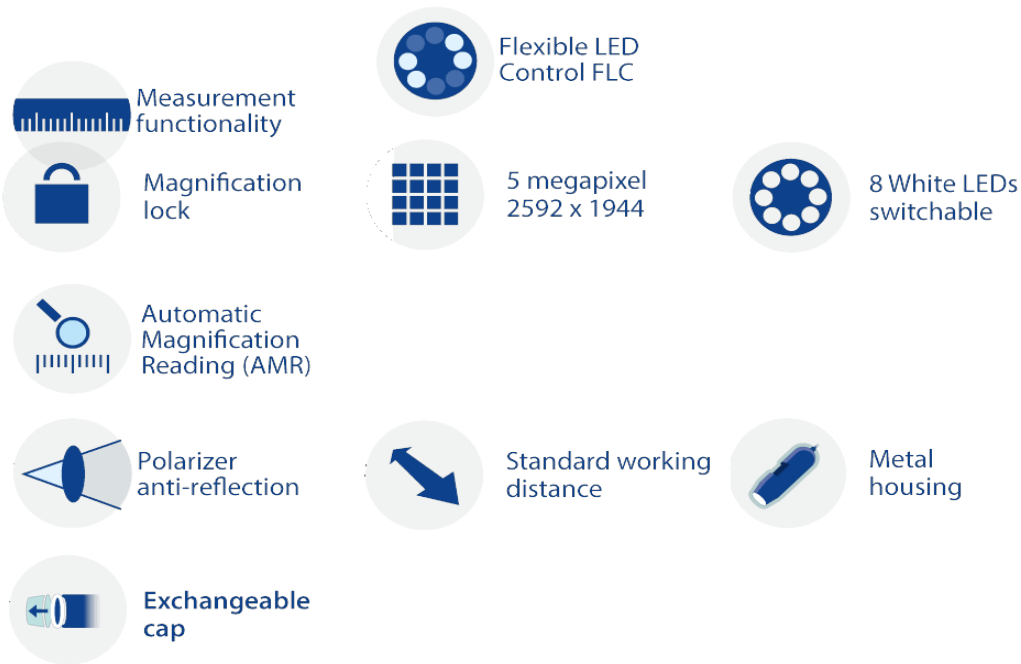

The Dino-Lite AM73515MZT is a USB 3.0 model with a magnification range of 10x - 220x and a standard working distance of up to 6 cm. This model is member of the Dino-Lite Edge family with the latest, cutting-edge optics, an unmatched 5 megapixel sensor and several special features such as Automatic Magnification Reading (AMR), Flexible LED Control (FLC), built-in polarization filter, full metal body, extensive measurement functions, calibration, exchangeable front caps and more.

Working distance/field of view/depth of field

MAGNIFICATION RATE	WORKING DISTANCE*	FIELD OF VIEW(X)	FIELD OF VIEW(Y)	DEPTH OF FIELD
10	142.4	37.8	28.3	19.7
20	59.5	19.5	14.6	5.4
30	32.9	13.0	9.7	2.6
40	20.2	9.8	7.3	1.6
50	13.2	7.8	5.8	1.1
60	9.0	6.5	4.8	0.8
70	6.4	5.6	4.2	0.6
80	4.8	4.9	3.6	0.5
90	3.9	4.3	3.2	0.42
100	3.5	3.9	2.9	0.36
110	3.3	3.6	2.7	0.31
120	3.5	3.3	2.4	0.27
130	3.8	3.0	2.2	0.24
140	4.3	2.8	2.1	0.22
150	4.9	2.6	1.9	0.2

160	5.6	2.4	1.8	0.18
170	6.4	2.3	1.7	0.17
180	7.3	2.2	1.6	0.15
190	8.2	2.1	1.5	0.14
200	9.2	2.0	1.5	0.13
210	10.2	1.9	1.4	0.13
220	11.3	1.8	1.3	0.12
Listed values may differ slightly	*Without front cap			Unit = mm

Specification

Lighting	
Light/ LED type	White
Number of LEDs	8
LED on/off switchable:	Yes
Infrared filter	IR cut-filter >650 nm
Diffuser available	Yes (N3C-D included)
Emission filter	No
Polarizer	Yes, linear
Optics	
Magnification	10-220x
Macro zoom	No
Working distance	Standard
Lens type	Glass with anti-reflection coating
Sensor	
Sensor type	CMOS
Resolution	5 Megapixel (2592x1944)
Maximum frame rate	45fps (max 20fps video recording)
Compatibility	
Interface	USB 3.0, Cable included**
Operating system	Windows7, 8, 10 & 11, MacOS 10.9 and up
Included software	DinoCapture 2.0 (Windows), DinoXcope (Mac OS)
Supported image formats (Windows)	BMP, GIF, PNG, JPG, TIF, RAS, PNM, TGA, PCX, MNG, WBMP, JP2, JPC, PGX
Supported video formats (Windows)	WMV, FLV, SWF
Supported image formats (MacOS)	JPEG, PNG
Supported video formats (MacOS)	MOV (max 1.3MP)

Imaging standards	DirectShow, UVC
Wifi	No
Housing	
Housing material	Metal housing
Magnification lock	Yes
Dimensions	11.9cm (L) x 3.3cm (H)
Weight	110g
Cable length	1.8m
Features	
Special feature	Automatic Magnification Reading (AMR) Flexible LED Control (FLC)
Measurement	Yes
Calibration	Yes
Microtouch sensor	Yes
ESD safe	Yes
Information	
Package contents	Microscope, carry pouch, software CD, user manual, quick guide, calibration target, front cover N3C-O- Open cap, N3C-C- Closed cap, N3C-D- Diffuser cap, N3C-E- Extension cap, N3C-L- Long cap, N3C-S- Side light cap
Warranty information	2 years European warranty
Regulatory approval	CE, FCC, ROHS
Price range	€1100,00 - €1300,00
Note	**Dino-Lite USB 3.0 models use a custom USB-A and USB-C cable that is specifically designed for Dino-Lite

Product Gallery